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
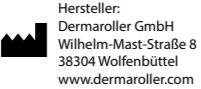
HC902 09.2024

HC902 Microneedling Roller The Original - 100% Made in Germany

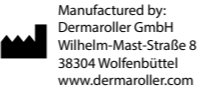Manufactured by:
Dermaroller GmbH
Wilhelm-Mast-Straße 8
38304 Wolfenbüttel
www.dermaroller.com

Thank you for choosing this high-quality product.

Using it two to three times a week will give your skin a fresher look and glow and it will visibly improve your complexion.

Dermaroller GmbH is a German manufacturer of medical and aesthetic devices and cosmetics. We have been manufacturing microneedling devices at our own production site since 1999.

That's why we call our microneedling roller "The Original".

Mode of Operation and Design:

Use your HC902 Microneedling Roller on smaller skin areas like face, neck and décolleté. The roller head has 162 tiny stainless-steel microneedles with an effective length of 0,2 mm.

When rolling across the skin the microneedles penetrate the corneal layer of the skin and remove loose skin cells and deposits. The result will be a fresher and more radiant looking skin.

Make the treatment with the HC902 part of your regular skin care program. This quality product is re-usable many times with proper care - see care instruction below.

WARNING: Keep out of reach of children.
CAUTION: Sharp needles - handle with care. Do not use this product near eyes to avoid injury.

Instruction for use

Please read this instruction carefully. It contains important information on the proper and safe handling of your HC902 device. Keep these instructions with the product as a reference.

Safety Notes:

- If used according to this instruction, harmful side effects are unknown.
- The device is to be used on adults > 18 years.
- Do not use on mucous membranes and lips.
- For hygienic reasons, the product must not be shared between individuals.
- The skin to be treated must be healthy. In case of an active skin disease or disorder (acne, eczema, etc.) seek professional advice first.

- Check your device for any damage before each use. Do not use a damaged device. Dispose of it safely in that case.
- Excessive pressure may cause temporary painful redness of the skin. It may also damage the needles.
- Do not soak the roller head in alcohol or solutions containing chlorine ions as this may damage the needles.
- Do not use harsh detergents, solvents or household cleaners. Do not sterilize with heat or steam. Do not clean in an ultrasonic bath. High temperatures > 185 °F / 85 °C affect the sharpness of the needles and can damage the plastic parts.

Use instruction:

- Clean the skin area to be treated e.g. with a skin cleanser and let it dry.
- Roll the device with gentle pressure 2 to 3 times in each direction (vertical, horizontal, diagonal left and right) across the skin in a star shaped manner to achieve a uniform application pattern. You may stretch the skin gently between your fingers when rolling across. Larger skin areas should be divided into smaller areas to allow for uniform results.
- Let the skin recover for 10 -15 minutes.
- Afterwards apply a high-quality moisturizing cream or lotion that contains no skin-irritating ingredients.

Tip: Especially sensitive skin may react with redness and scaling after the first few applications. This is a normal, temporary reaction. Start with only a single application per week until your skin gets used to it. Afterwards you can increase frequency of application.

Care Instruction - Careful treatment and cleaning of the product can extend product life:

- Rinse the device under plenty of hot, but not boiling tap water and shake it dry. Avoid contact with hard surfaces.
- We advise to disinfect the needles with a suitable disinfectant. Spray the roller head all around with 4 puffs of a mild disinfectant/bactericidal solution that is approved for use on skin. Let it dry for a few minutes.
- Afterwards, put your HC902 back in the provided packaging (Jewel Case or Blister) and carefully store it in a safe place.

Product life:

All sharp steels are subject to wear, even when "only" applied to skin. For the HC902 a product lifetime of at least 50 applications can be expected when used and cleaned as recommended above. Microneedling of larger skin areas may wear off the needles faster and shorten product life.

HC902 Microneedling Roller is a consumer product that meets the relevant safety requirements.

This product meets CPSC safety requirements.

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